

Siemens Healthcare GmbH, HC AT IR D, Siemensstr. 1, 91301 Forchheim

Name	Daniel Seebacher
Department	HC AT IR D
Telephone	+49 (9191) 18-6427
Fax	+49 (9191) 18-9920
Mobile	+49 (173) 4310033
E-mail	daniel.seebacher@siemens.com
Date	January 24, 2017

**To all users of Artis systems with software versions
VC14, VC21, VD10, VD11 and VE10**

Important customer safety notice regarding corrective field action:

AX058/16/S

**Information regarding corrective action for Artis systems with software versions VC14, VC21, VD10,
VD11 and VE10**

Dear Customer,

This letter is to inform you about an existing system function that is not fully described in the operating instructions. The function concerned may allow you to complete a treatment or diagnosis despite a defective focus. For this reason, we consider the incompleteness of the operating instructions to be safety-relevant.

What is the underlying issue requiring this corrective action and when does the issue occur?

For Artis systems with the software versions specified above, a semi-automatic focus switch is implemented. Interacting with a special type of defect due to wear in the foci "small" and "micro" (if present), this performs an automatic focus switch when the footswitch is pressed several times. This special system function is not fully described in the operating instructions currently available. If a focus (small/micro) develops a defect, an alternative focus can be selected by pressing the footswitch several times. This enables radiation to be generated once again so that the clinical treatment or diagnosis can be continued.

The attached supplement to the operating instructions describes the corresponding "system messages" and explains how to use the function properly.

Siemens Healthcare GmbH
Management: Bernhard Montag, Chairman;
Thomas Rathmann, Michael Reitermann

Siemensstr. 1
91301 Forchheim
Germany

Tel.: +49 (9191) 180
siemens.com/healthcare

Chairman of the Supervisory Board: Siegfried Russwurm
Registered offices: Munich, Germany; Commercial registry: Munich, HRB 213821
WEEE Reg. No. DE 64872105

What action can you take?

Please bring the attached supplement to the operating instructions to the attention of your staff responsible for operating the Artis systems, and ensure that the information is available in the event of a fault.

Thank you for your cooperation in connection with this customer safety notice. Please immediately inform and instruct all employees in your organization who need to be aware of this problem. Please also forward this safety information to any other organizations that could be affected.


If the affected device has been sold and is therefore no longer in your possession, please forward this safety notice to the new owner. We would also ask you to inform us of the identity of the device's new owner where possible.

Best regards,

Siemens Healthcare GmbH
Advanced Therapies Business Area



Dr. Heinrich Kolem
President Advanced Therapies

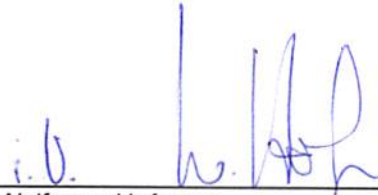

Wolfgang Hofmann
Safety Officer Medical Devices